

Accessing Dexcom Clarity from your computer and uploading your receiver

Dexcom clarity is a cloud-based depository for glucose data which allows you to look at your past glucose data, glucose patterns and trends.


1. Open web browser and go to <https://www.clarity.dexcom.eu/>
2. Click on “Dexcom Clarity for home users”
3. Put in your username and password and click “Login”
4. If you have not uploaded your receiver from your computer at home before, the following instructions will appear – follow the steps

We Don't Have Your Data Yet

Follow these steps to upload your data and see it displayed through Dexcom CLARITY.
If you've previously completed Steps 1 and 2 on this computer, you can jump to Step 3.

Step 1

Download Dexcom CLARITY Uploader

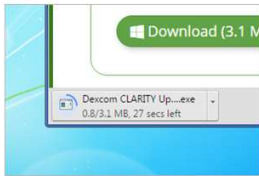


Download the Dexcom CLARITY Uploader, which will connect your CGM receiver with Dexcom CLARITY.

Download (3.7 MB)

Step 2


Open Downloaded File



When the file has downloaded, open it and follow the installation instructions.

Step 3

Connect Your CGM Receiver

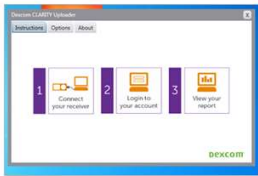

Plug your CGM receiver into your computer and follow the Dexcom CLARITY Uploader instructions.

5. For all future uploads, login to Dexcom clarity using the above instructions and just connect your receiver to your computer.
6. **Dexcom receivers only store about 30 days of CGM readings. It is very important to upload your Dexcom receiver every 3 weeks to make sure we don't lose data.**